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## INTRODUCTION

*Dear customer, thank you for choosing our product. We feel certain you will appreciate its quality and reliability as it was designed and manufactured with customer satisfaction in mind. These instructions for use comply with the European Standard EN 82079.*

### **WARNING! Instructions and warnings for a safe use**

**Before using this appliance, carefully read the instructions for use and in particular the safety notes and warnings, which must be complied with. Keep this manual together with its illustrative guide for future consultation. Should you transfer the appliance to another user, make sure to hand over this documentation as well.**

**NOTE:** should any part of this manual be difficult to understand or any doubt arise, contact the company at the address indicated on the last page before using the product.

**NOTE:** partially sighted persons can consult the digital version of these instructions for use on the website [www.tenactagroup.com](http://www.tenactagroup.com)

## INDEX

Safety notes	page 33
Symbols	page 37
Product use contraindications	page 38
Other precautions for safe use	page 42
Possible side effects	page 45
Description of the appliance and its accessories	page 48
Treatment areas	page 49
Hair removal program	page 50
Information before use	page 51
Use	page 55
Maintenance	page 58
Troubleshooting	page 58
Disposal	page 59
Assistance and warranty	page 60
Skin tone range	I
Product technical data	II
Adaptor technical data	II
Illustrative guide	III

## SAFETY NOTES

- After removing the appliance from its packaging, verify its integrity based on the drawing and the presence of any damage caused by transport. If in doubt, do not use the appliance and contact an authorised service centre.

**CAUTION!** Remove any communication material such as labels, tags etc. from the appliance before use.

**WARNING!** Risk of suffocation.

Children shall not play with the packaging. Keep the plastic bag out of reach of children.

- This appliance must be used only for the purpose for which it was designed, namely as an epilator for home use. Any other use is considered not compliant with the intended use and therefore dangerous.
- Before connecting the appliance to the mains supply, make sure that the voltage and frequency indicated on the technical data of the appliance correspond to those of the mains supply. These technical data are present on the appliance or on the power supply unit, if any (see illustrative guide).
- This appliance can work automatically at a frequency of 50 Hz or 60 Hz.
- The appliance is only to be used with the power supply unit provided.
- This appliance can be used by persons aged from 16 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or

instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not use the appliance. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children.

- DO NOT use the appliance with wet hands, damp feet or barefoot.
- DO NOT pull the supply cord or the appliance itself to remove the plug from the socket.
- DO NOT pull or lift the appliance by the supply cord.
- DO NOT expose the appliance to humidity or atmospheric agents (rain, sun, etc.) or extreme temperatures.
- The appliance must not be left unattended while it is connected to the mains supply.
- Always disconnect the appliance from the mains supply before cleaning or performing maintenance and when the appliance is not in use.
- If the appliance is faulty or malfunctions, switch it off and do not tamper with it. Contact an authorised technical service centre for any repairs.
- If the supply cord is damaged it must be replaced by an authorised technical service centre, so as to prevent all risks.

- When the appliance is used in a bathroom, unplug it after use since the proximity of water presents a hazard even when the appliance is switched off.
- For additional protection, the installation of a residual current device (RCD) having a rated residual operating current not exceeding 30 mA (ideal 10 mA) is advisable in the electrical circuit supplying the bathroom. Ask your installer for advice.

**WARNING! DO NOT** use this appliance near bathtubs, showers, basins or other vessels containing water.

- This appliance shall not come into contact with water, other liquids, spray, steam. For any cleaning and maintenance operation refer to the instructions in the relevant chapter.
- DO NOT use this appliance on animals.
- This appliance is for personal use only: DO NOT share the use of the appliance.
- For the features of the appliance, refer to the external packaging.
- Use the appliance only with the supplied accessories that are an essential part of the appliance.
- This appliance is NOT suitable for use in the presence of flammable

anaesthetic mixtures with air or with Oxygen or Nitrous Oxide.

- DO NOT obstruct the ventilation filters.
- DO NOT use the appliance if the skin colour sensor is cracked or broken.
- DO NOT use the appliance if the external casing is cracked or broken
- DO NOT use the appliance if the lamp is damaged or if the treatment surface is cracked, chipped or missing.

## SYMBOLS

	<i>Warning</i>		<i>Double insulation</i>
	<i>Prohibited</i>		<i>Only for indoor use</i>
	<i>Not suitable for use in the bath or shower</i>		<i>Away from the eyes</i>
	<i>Stand-by</i>		<i>Device ready for pulse emission</i>
	<i>Power supply connector polarity</i>		<i>No tanned skin</i>

## PRODUCT USE CONTRAINDICATIONS

**WARNING!** Safety notes. Read before using the appliance.

**This appliance is NOT suitable for use by everyone. Read and take all information contained in the following paragraph into account before using the product. For more information and personalised advice visit [www.bellissima.com](http://www.bellissima.com).**

If you are unsure whether you can use the appliance or not, ask your doctor or dermatologist!

**⊘ DO NOT USE THE APPLIANCE IF ONE OF THE FOLLO WING CONDITIONS APPLIES:**

- During pregnancy or nursing
- After exposure to the sunlight or artificial tanning in the past 28 days
- If you have tattoos or permanent makeup in the treatment area
- If you have dark brown or black spots, such as large freckles, birth marks, moles or warts in the treatment area
- If you have active eczema, psoriasis, lesions, open wounds or infections (such as herpes outbreaks) in the treatment area. Wait for the area to be healed before using the appliance
- If you have a history of keloidal scar formation
- If you are aware of being photosensitive
- If you take medicines which make skin more sensitive to light,

including non steroidal anti-inflammatories (like aspirin, ibuprofen, paracetamol), tetracyclines, phenothiazines, thiazide, diuretics, sulfonylurea, sulfonamides, DTIC, fluorouracil, vinblastine, griseofulvin, alfa hydroxy acids (AHA), beta hydroxy acids (BHA), Retin-A®, Accutane® and / or topical rethinoids or azelaic acid.

- If you have abnormal skin conditions, like albinism or vitiligo
- If you have varicose veins in the treatment areas
- If you are being or have recently been treated with alfa hydroxy acids (AHA), beta hydroxy acids (BHA), Retin-A®, topical rethinoids or azelaic acid
- If you have abnormal skin conditions caused by diabetes or other systemic or metabolic disorders
- If you have been treated with Accutane® (isotretinoin) for the past 6 months
- If you have been on steroids during the last 3 months
- If you have a history of herpes outbreaks in the treatment area, unless you have consulted your doctor and have been treated before using the appliance
- If you suffer from epilepsy
- If you carry an active system, such as a pacemaker, an incontinence aid, an insulin pump, etc.
- If you suffer from diseases related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc.
- If you have a history of skin cancer or areas of potential skin malignancies
- If you have undergone radiation therapy or chemotherapy within the past 3 months

- **Any other problem or condition which in your physician's opinion would make treatment unsafe for you.**

**WARNING! The appliance may be less effective in removing hair from women with hormonal disorders since they might cause excessive hair growth, hirsutism or increased hair growth on the face and neck. If in doubt, consult your doctor.**

#### THE COLOUR OF YOUR SKIN

The appliance removes unwanted hair and selectively affects the hair pigment. Various amounts of pigment (melanin) are also found in the surrounding skin tissue. The amount of pigment in the skin, which affects the skin tone, determines the risk posed to the skin when any type of pulsed light epilation is used.

**WARNING! Before using the appliance, compare your complexion with the skin tone range provided in this manual and/or on the package of the product the appliance.**

#### NATURALLY DARK SKIN

**WARNING! DO NOT use the appliance on naturally dark or tanned skin.**

Using the appliance on dark skin can cause adverse effects such as burns, blisters and changes in the colour of the skin (hyper or hypo-pigmentation).

#### SKIN COLOUR SENSOR

The appliance has a skin colour sensor (1B) that measures the tone of the skin in the area that is to be treated at the beginning of each session and occasionally during the session. The skin colour sensor (1B) ensures that pulses are only emitted on suitable skin tones.

#### TANNED SKIN

**DO NOT use the appliance with energy levels from 2 to 5 on tanned skin or after recent exposure to the sun and/or tanning beds. Do not expose the treated areas to the sun. Such use can cause serious burns, wounds, blisters and changes in the colour of the skin (hyper or hypo-pigmentation).**

Tanned skin, especially after sun exposure, contains large amounts of melanin.

This applies to all types of skin and complexions, including those that do not seem to tan quickly. The presence of large amounts of melanin, expose the skin to greater risk when light epilation methods are used.

**If the appliance is going to be used at energy levels between 2 and 5, avoid exposure to the sun and/or tanning beds for 4 weeks**

**before the treatment. The appliance can be used after exposure to the sun only if the energy level is set to the lowest level (level 1).**

Exposure to the sun means constant unprotected direct exposure to sunlight for more than 15 minutes or constant unprotected indirect exposure to sunlight for more than 1 hour.

Always protect recently treated skin when exposed to the sun, making sure to carefully spread sun cream SPF 30 or higher for 2 weeks after each hair removal session.

## OTHER PRECAUTIONS FOR A SAFE USE

### THE EYES

**DO NOT use the appliance around or near the eyes.**

The appliance can be used to remove hair in a delimited part of the face. Special care must be taken to avoid the eyes when treating this area.

### BEFORE THE TREATMENT

- **DO NOT look directly at the light emitted from the treatment surface.**
- **DO NOT point the treatment surface whilst emitting pulsed light towards an empty space.**

**WARNING! Keep away from the eyes.**

- **DO NOT use the appliance on any part of the body where you may wish hair to grow subsequently.**
- **DO NOT USE THE APPLIANCE ON NIPPLES OR GENITAL ORGANS.**

### THE POWER LEVEL

Select the power level CAREFULLY!

The power level refers to the intensity of the pulse of light that is projected on the skin and ranges from the lowest to the highest. The power level LEDs (5B) indicate the set power level.

When the power level increases, the results improve, however, the risk of side effects also increases (refer to “Possible side effects”).

- **Always start the treatment at the lowest energy level.**
- **Always shave the area that is to be treated and make sure the skin is clean and dry before using the appliance.**

The best results are achieved 3 days after shaving, when the hairs are 1-2 mm long.

Only if you feel little or no discomfort during and after use with the power level set to the lowest can you increase the power by one level, and so forth until you reach the maximum level suitable for your skin type. For detailed instructions regarding the power level, refer to “Power Level” in “Information before use”.

- On your face, if there are few and thin dark hairs to be treated, it is not necessary to shave the area, only trim them shorter with scissors. If the hairs are very short (approximately 1-2 mm) and not too thick, there is no need to shave nor trim them.

- **Cover any birthmark and tattoos before using the appliance.**
- **Cover dark brown or black spots, such as large freckles, a birthmark, moles or warts before using the appliance.**
- **⊘ DO NOT use flammable liquids, such as alcohol, products containing alcohol (including perfume, sanitary products, etc.) or acetone to clean the skin before using the appliance.**

#### THE PATCH TEST

It is recommended to perform a patch test before starting treatment, in order to identify the maximum tolerance of your skin to light pulses. Refer to “Information before use” for instructions regarding the patch test.

#### DURING TREATMENT

**WARNING! AVOID treating the same skin area more than once or overlapping skin areas.**

**WARNING! DO NOT overlap pulses.**

- **If the skin blisters or burns, STOP USE IMMEDIATELY.**
- **Ensure that the treatment surface is pointed towards the skin, and is completely in contact with it.**
- **Remove the treatment surface from the skin if the skin or the treatment surface become excessively hot.**

#### AFTER TREATMENT

- After treatment it is possible to apply an alcohol and perfume-free, gentle or soothing moisturising cream or lotion to the area.
- **After treatment at energy levels between 2 and 5, DO NOT intentionally expose the skin to sunlight and/or solarium lamps.**
- **Apply sunscreen protection SPF 30 or higher to skin treated at any energy level when not intentionally exposed to sunlight for 2 weeks after each hair removal session.**

**WARNING! The use of the appliance may cause temporary changes in pigmentation. (See “Possible side effects” section).**

#### POSSIBLE SIDE EFFECTS

Side effects and complications associated with the appliance are rare when used according to the instructions.

However, any cosmetic procedure, including those intended for home use, pose a potential risk. It is therefore important to understand and acknowledge the risks and complications that may occur when using pulsed light epilation systems designed for home use.

**In case of any side effects or complications, contact your doctor or dermatologist.**

Negative Event	Probability of the Negative Event	Negative Effect	Probability of the Negative Effect
Accumulated multiple pulses or overlapping on the same section of skin	Minimum	Minimal Discomfort to the Skin	Minimum
		Redness on the Skin	Minimum
		Greater Skin Sensitivity	Minimum
		Wounds and burns	Rare
		Scars	Rare
		Changes in Pigmentation	Rare
		Excessive Swelling and Redness	Rare
		Infection	Negligible
		Bruising	Negligible

### Minimal Discomfort to the Skin

Despite home use pulsed light epilation is generally very well tolerated, most users experience mild discomfort during use, which is usually described as a slight feeling of warmth or tingling in the treated sections. This feeling can be felt during the treatment or a few seconds after.

Any other discomfort is not normal and treatment must be stopped because your skin does not tolerate the hair removal treatment or the power level setting is too high.

### Redness on the Skin

The skin may become red immediately after using the appliance or within the subsequent 24 hours. The redness usually disappears within 24 hours. Consult your doctor if the redness does not disappear within 2 or 3 days.

### Greater Skin Sensitivity

The treated skin becomes more sensitive and may become dry or flaky.

### Wounds and Burns on the Skin

Very rarely, burns or wounds may appear after the treatment. Burns or wounds can take several weeks to heal and very rarely, may leave a visible permanent scar.

### Scars

Even though it is very rare, permanent scarring may occur. Generally, when scars are formed, they appear flat and white (hypotrophic).

However, the scar may also be big and red (hypertrophic) or big and beyond the wound itself (keloid). Subsequent aesthetic treatment may be required to improve the appearance of the scar.

### Changes in Pigmentation

In some rare cases, there lies the risk of temporary hyperpigmentation (increased pigment or brown discolouration) or hypopigmentation (whitening) in the skin surrounding the hair. The risk of changes in skin pigmentation is higher in persons who have a dark skin tone. Generally, discolouration or changes in the skin pigment are temporary, and permanent hyperpigmentation or hypopigmentation is rare.

### Excessive Swelling and Redness

In some rare cases, the skin may become very red and swollen. This is more common in sensitive parts of the body. The redness and swelling should lessen in the following 2 to 7 days and ice must be applied often. The area can be cleaned gently, however, avoid exposure to sun.

## Infection

Skin infection is very rare, however, it may still occur following a burn or wound caused by the appliance.

## Bruising

Very rarely, blue-purple bruising may be caused after using the appliance and may last for 5 to 10 days. As the bruise fades, rust-brown skin discolouration may occur (hyperpigmentation), which may be permanent.

Please consult your doctor or dermatologist in the case of side effects or complications.

## DESCRIPTION OF THE APPLIANCE AND ITS ACCESSORIES

Home use pulsed light appliance, based on HPL technology (Home Pulsed Light) to permanently reduce unwanted hair. Permanent hair reduction means long-term reduction in the amount of hairs that grow after a course of treatment.

The appliance emits a light pulse that is absorbed by the pigment in the hair. The light beam is conveyed to the hair follicle, thereby generating heat that can deactivate it. The appliance is also equipped with two electrodes that generate galvanic energy to boost the action of the light pulse.


Refer to the Illustrative guide (III) to check the content of the package.

1. Lamp
  - 1A. Treatment surface
  - 1B. Skin colour sensor
  - 1C Galvanic energy contact electrode
2. Pulse activation button
3. Power supply unit connection
4. Ventilation filters

5. Control panel
  - 5A. Power level setting and on/off button
  - 5B. Power level LEDs
  - 5C. LEDs to check the skin colour and indicate the device is ready.
6. Product technical data
7. Galvanic energy base electrode
8. Plug
  - 8A. Plug A, 8B. Plug B, 8C. Plug C
9. Plug technical data
  - A. Skin tone range

For the appliance features, refer to the external packaging.  
Use only the supplied accessories.

## TREATMENT AREAS

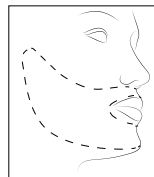The appliance can be used to remove unwanted body hair in the following areas:  
underarms, bikini line, arms, legs, back and chest.

Women can also use it in the part indicated by the white dashed line.

It is important to only treat the facial areas where there is visible thick and dark superfluous hair. Areas with thin fair hair, often on and under the chin, must not be treated.

## HAIR REMOVAL PROGRAM

A typical cycle of hair growth can take 18 to 24 months. During this time, more sessions may be required to achieve permanent hair removal.

Hair removal efficacy varies from person to person and in different body parts, according to the hair colour and level of tolerance of one's own skin.

Recommended treatment plan in a complete hair growth cycle:

- From session 1 to session 4: sessions must be held two weeks apart from each other.
- From session 5 to session 7: sessions must be held four weeks apart from each other.
- The subsequent maintenance sessions can be performed from time to time, when and if necessary, until long-lasting results are achieved.

**ⓘ NOTE: The hair removal results will not be enhanced if the same part of the body is treated more often than recommended (every two weeks).**

Dormant hair follicles can be reactivated due to hormonal and physiological changes. Therefore, occasional hair removal sessions may be required.

The individual reaction depends on the type of hair and biological factors that may affect the growth cycle of hair. Some people can react quicker or slower with regards to the recommended hair removal plan.

In any case, a long-lasting result cannot be expected with only one or two sessions.

**The appliance is not effective on white hair. The appliance may be less effective on natural grey, light blonde and red hair.**

## INFORMATION BEFORE USE

### POWER LEVEL

The appliance can be set to 5 power levels; 1 (lowest) to 5 (highest), indicated by 5 power level LEDs (5B).

### SETTING THE POWER LEVEL

Press the power level setting and on/off button. Each time it is pressed, the power level increases by 1 up to 5, after which it returns to level 1.

Perform this action until the desired level is set.

The green power level LED (5B) will light up in line with the set power level.

**WARNING! Always start a hair removal session by setting the energy level detected by the patch test.**

### SKIN COLOUR SENSOR (1B)

Using the pulsed light epilation system on dark skin can cause adverse effects such as burns, blisters and changes in the colour of the skin (hyper or hypo-pigmentation).

The appliance has a skin colour sensor (1B) that measures the tone of the skin in the area that is to be treated at the beginning of each session and occasionally during the session.

If the skin colour sensor (1B) detects a complexion that is too dark for the treatment, the appliance stops emitting the light pulses automatically and the skin colour/device-is-ready LED (5C) starts to flash rapidly.

If this occurs, try to use the appliance on another part of the body or contact the IMETEC service centre.

## THE TWO APPLICATION METHODS

The device can be applied in two ways: “by Pulse” and “Sliding”.

### By Pulse Method

The by pulse method is optimal when using higher power levels as it allows for better control of the skin areas being covered. Each time the pulse activation button is pressed (2), a light pulse is emitted.

A precise technique must be followed when using the by pulse method to make sure the entire treatment area is covered and dangerous overlapped pulses are prevented:

- Imagine the treatment area is divided in rows.
- Emit a light pulse proceeding in order, starting from one end of each row and continuing to the other end (from the bottom upwards or right to left, depending on the treatment area).

The area covered by the treatment surface (1A) leaves temporary marks on the skin due to the light pressure applied by the handpiece.

These marks remain visible for a few moments and are useful to establish the subsequent pulse position.

### Sliding Method

The Sliding method is optimal when using low power levels and it allows for faster and more efficient sessions. When using the Sliding method, keep the pulse activation button (2) pressed and slide the device in a slow gliding motion while remaining in contact with the skin, along the area to be treated.

At lower power levels (1 and 2), the device takes only 1 to 1.5 seconds to recharge between one pulse and another, thereby allowing an uninterrupted sequence of flashes.

**WARNING! DO NOT overlap impulses.**

**DO NOT treat the same area of skin more than once per hair removal session!**

**If an excessive heat is felt, skin itches or it blisters or burns, STOP USE IMMEDIATELY.**

### HOW TO PERFORM THE PATCH TEST

Since the skin tone and its sensitivity vary from one body part to another, a patch test should be performed on each part of the body that is to be treated.

- First check the tone of your skin by comparing it to the table on the package and in the instruction manual.
- Identify a small area that is to be treated.
- Shave the area that is to be treated.
- After setting the lowest power level, apply the handpiece to the area. Grip the handpiece correctly making sure that your hand is touching the galvanic energy base electrode (7). If you do not feel any discomfort, deliver the first pulse. If no discomfort is felt after the pulse, increase the power to the next level, move the handpiece to the adjacent skin area and emit another pulse. Increase the power level on each subsequent area of skin until the maximum tolerable level for the skin type. A feeling of heat can be felt on the skin during the treatment, and this is normal. If discomfort is felt, the maximum tolerable power level is the previous one.

- Check the area after 1 hour: if there is no reddening and/or discomfort, it is possible to start treatment in the tested area at the last energy level setting. If redness persists, the maximum tolerable power level for your skin type is the previous one. Repeat the patch test on every part of your body that you wish to treat.

Even though it is extremely rare, if hypersensitivity should persist for over 24 hours, it is recommended to consult a doctor.

**WARNING! The power output directly affects the result, however, due to the principle of heat of the treatment, the maximum limit tolerated by one's skin must not be exceeded. This limit is normally high for fair skin and low for dark skin because melanin, that darkens the skin, increases the production of heat. DO NOT exceed the maximum limit that your skin tolerates with no discomfort whatsoever.**

#### PREPARING THE AREAS THAT ARE TO BE TREATED

- The skin must be shaved, clean, dry, with no make-up, anti-perspirant or deodorant.
- On your face, if there are few and thin dark hairs to be treated, it is not necessary to shave the area, only trim them shorter with scissors. If the hairs are very short (approximately 1-2 mm) and not too thick, there is no need to shave nor trim them.

- Moles, birthmarks, dark spots, tattoos and lips must be covered with a white pencil or white patches.


**WARNING! Carefully read all warnings and safety precautions, possible adverse reactions, treatment areas and epilation programme.**

## USE

#### PREPARING THE APPLIANCE

1. Remove the appliance and the other components from the package.
2. If the appliance is stored in a very cold place, wait for it to reach at least 15 °C.
3. If the plug has not been inserted into the power supply unit (8), hold the power supply unit, pull the "Open" key downwards and insert the plug (8A, 8B or 8C) which is suitable for the socket of the electrical system, into the relative slot.
4. Insert the cable of the plug (8) into the power supply connection (3).
5. Insert the plug (8) into a socket: all the LEDs on the control panel (5) will go on and off repeatedly in a circular sequence.

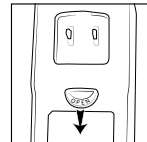**WARNING! DO NOT use this appliance if there are visible signs of damage.**

**TREATMENT**

6. Press the power level setting and on/off button (5A): the green power level 1 LED (5B) lights up and the ventilation system is activated (you will hear a sound similar to that of a hair dryer).
7. About 1 second after having pressed the power level setting and on/off button (5A), the skin colour/device-is-ready LED (5C) lights up: the device is automatically set to the lowest power level.
8. Set the energy level that you wish to use for the treatment. To set the energy level, follow the directions provided in the 'Information before use' chapter.
9. Apply the treatment surface (1A) firmly on the skin. Only when the treatment surface (1A) and the galvanic energy base electrode are entirely in contact with the skin, the skin colour/device-is-ready LED (5C) starts to flash slowly. Be sure to touch the galvanic energy base electrode (7) located on the back of the appliance, with your hand, to activate the galvanic energy circuit.
10. Press the pulse activation button (2). The appliance will immediately check your skin colour: if the skin colour is fair enough for safe application the appliance will activate the galvanic energy circuit and you may experience a slight tingling sensation. This is a normal sensation: galvanic energy helps open the pores, improving the action of the light pulse. Immediately afterwards, the appliance delivers the light pulse onto your skin, accompanied by a sound, and the skin colour/device ready indicator light (5C) switches off. A slight tingling sensation may be felt. Stop using the appliance immediately if excessive heat or tingling is felt or blisters or small burns appear on the skin.  
 Note: if the appliance does not emit a light pulse and the skin colour/device-is-ready LED (15) flashes rapidly, the skin tone is too dark for safe application. Try to use the appliance on another part of the body.
11. Remove the treatment surface (1A) from the treated area.

12. A few seconds (depending on the power level set) after the previous pulse is emitted, the skin colour/device-is-ready LED (5C) goes on again.
13. Select the application mode to continue the treatment.

**By Pulse Method**

- 14A. Move the treatment surface (1A) to another section. Use the pressure marks left previously to place the treatment surface on the adjacent section, avoiding overlapped pulses (See "Covering the entire treatment area").
- 15A. Repeat the process, starting from point 8, until the entire area to be treated is covered.

**Sliding Method**

- 14B. Keep the pulse activation button (2) pressed and slide the device in a slow gliding motion while remaining in contact with the skin, along the area to be treated. The recharging time between pulses allows for an uninterrupted sequence of flashes. Make sure the entire area you intend to treat is covered while avoiding to go over the same area a number of times during the same session.

**AFTER TREATMENT**

16. Press the power level setting and on/off button (5A): all the LEDs will flash sequentially.
17. Remember the power level used for the treatment session as this will be reset when the appliance is next switched on.
18. Disconnect the adaptor cable (8) from the socket.

**Caring for the skin after the hair removal session.**

See "Other precautions for safe use".

## MAINTENANCE

**WARNING! Always disconnect the appliance from the mains supply before cleaning or performing maintenance.**

### CLEANING AND STORING

It is recommended to clean the appliance and the treatment surface (1A) after each hair removal session:

- Disconnect the appliance from the power supply (8) before cleaning it.
- Clean the surface of the handpiece and the treatment surface (1A) gently, with a dry, clean cloth.
- Store the appliance in a cool and dry place, away from heat and water.

## TROUBLESHOOTING

<b>The appliance does not go on</b>	<ol style="list-style-type: none"> <li>1. Make sure the plug (8) is inserted correctly into the power supply connection (3).</li> <li>2. Make sure the plug (8) is inserted into a wall socket.</li> </ol>
<b>When I press the pulse activation button (2) the galvanic pulse is not delivered</b>	<ol style="list-style-type: none"> <li>1. Make sure that the hand holding the appliance is touching the galvanic energy base electrode (7). To activate the galvanic energy it is necessary to press the pulse activation button (2) as your hand touches the galvanic energy base electrode (7) and the skin on the area to be treated is also touching the galvanic energy contact electrode (1C) at the same time.</li> </ol>

**No light is emitted when the pulse activation button (2) is pressed**

1. Make sure the skin colour/device-is-ready LED (5C) flashes, which indicates good contact with the skin, and the surface treatment (1A) is evenly and firmly pressed on the skin. OR
2. Verify the skin colour/device-is-ready LED (5C): if it flashes rapidly, the skin tone is too dark for safe application. Try to use the appliance on another part of the body or contact the IMETEC service centre. OR
3. Check the green power level LEDs (5B): a system error is present if one or more power level LEDs (5B) flash while the skin colour/device-is-ready orange LED (5C) remains off. If the problem persists, contact the IMETEC Customer Service Centre. OR
4. Verify whether the skin colour/device-is-ready orange LED (5C) and the green power level LEDs (5A) are flashing all at once. If this happens, disconnect the device from the power supply. The device has completed its operational life cycle and must be disposed of properly.

## DISPOSAL

The product is packaged in recyclable materials. Dispose of packaging in accordance with the environmental protection standards.

Pursuant to Art. 26 of Italian Legislative Decree no. 49 of 14 March 2014 "Implementation of Directive 2012/19/EU on waste electrical and electronic equipment (WEEE)", the crossed bin symbol on the appliance or its packaging indicates that the product must be disposed of separately from other waste at the end of its useful life.

The user must therefore take the appliance to an appropriate recycling centre for electrical and electronic equipment. Alternatively, instead of handling the disposal, the appliance can be handed over to the retailer when a new, equivalent appliance is purchased. Retailers of electronic products having a sales area of at least 400 m<sup>2</sup> can also be handed electronic products that are to be disposed of measuring

less than 25 cm, free of charge and with no obligation to buy. Appropriate separate collection for the subsequent procedure for the decommissioned appliance for recycling, which includes treatment and environmentally friendly disposal helps prevent adverse effects on the environment and health and promotes reusing and/or recycling materials that the appliance consists of.

## ASSISTANCE AND WARRANTY

The appliance is guaranteed for a period of 2 years from the date of purchase. Evidence shall be taken as the date indicated on the receipt/invoice (provided it is legible), unless the buyer can prove that the delivery occurred later. In the event of a defect in the product which existed prior to the date of delivery, appliance repair or replacement is guaranteed free of charge, unless one of the two solutions is disproportionate to the other. The buyer has the obligation to notify an authorised Service Centre of the non-conformity within two months from discovering the defect. The warranty does not cover any part that may be defective due to:

- a. Damage from transport or accidental drops,
- b. Incorrect installation or inadequate electrical system,
- c. Repairs or alterations made by unauthorised personnel,
- d. Poor or incorrect maintenance and cleaning,
- e. Product and/or product parts subject to wear and/or consumables,
- f. Failure to comply with the appliance's operating instructions negligent or careless use,

The list above is purely indicative and not exhaustive and in any case, this warranty shall not cover all **circumstances that cannot be attributed to manufacturing defects** of the appliance. In addition, the warranty does not cover any case of improper

use of the appliance and professional use. We disclaims any liability for any damage that may directly or indirectly be caused to persons, property and domestic animals as a result of failure to comply with all the instructions provided in this installation, use and maintenance booklet of the appliance.

This is without prejudice to any contractual warranty claims against the seller.

### Methods of support

Appliance repairs must be carried out by an authorised Service Centre. If the faulty appliance is under warranty, it must be sent to the Service Centre together with a fiscal document that can confirm the date of sale or delivery.

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© ELLISSIMA
Tenacta Group S.p.A. Via Piemonte, 5/11 24052 AZZANO S. PAOLO (BG) ITALY